



Distributor Benelux



**BLOOMEDICAL.**

PART OF ARSEUS MEDICAL GROUP

Neerloopweg 5A1, 4814 RS Breda, The Netherlands  
Tel. +31 76 524 53 53, [www.bloommedical.com](http://www.bloommedical.com)

Manufacturer

**EUFOTON**<sup>®</sup>  
medical lasers


**Eufoton Srl** Via Flavia 23/1 - 34148 Trieste - Italy  
Ph. +39 040 9899082 - Fax +39 040 9890555  
[www.eufoton.com](http://www.eufoton.com) - [info@eufoton.com](mailto:info@eufoton.com)  
[www.endolift.com](http://www.endolift.com)

**ISO  
9001**  
QUALITY  
ASSURANCE

**ISO  
13485**  
QUALITY  
ASSURANCE



THE BEST  
NON-SURGICAL  
LASER ASSISTED  
TREATMENT  
TO REDUCE  
CUTANEOUS  
LAXITY.

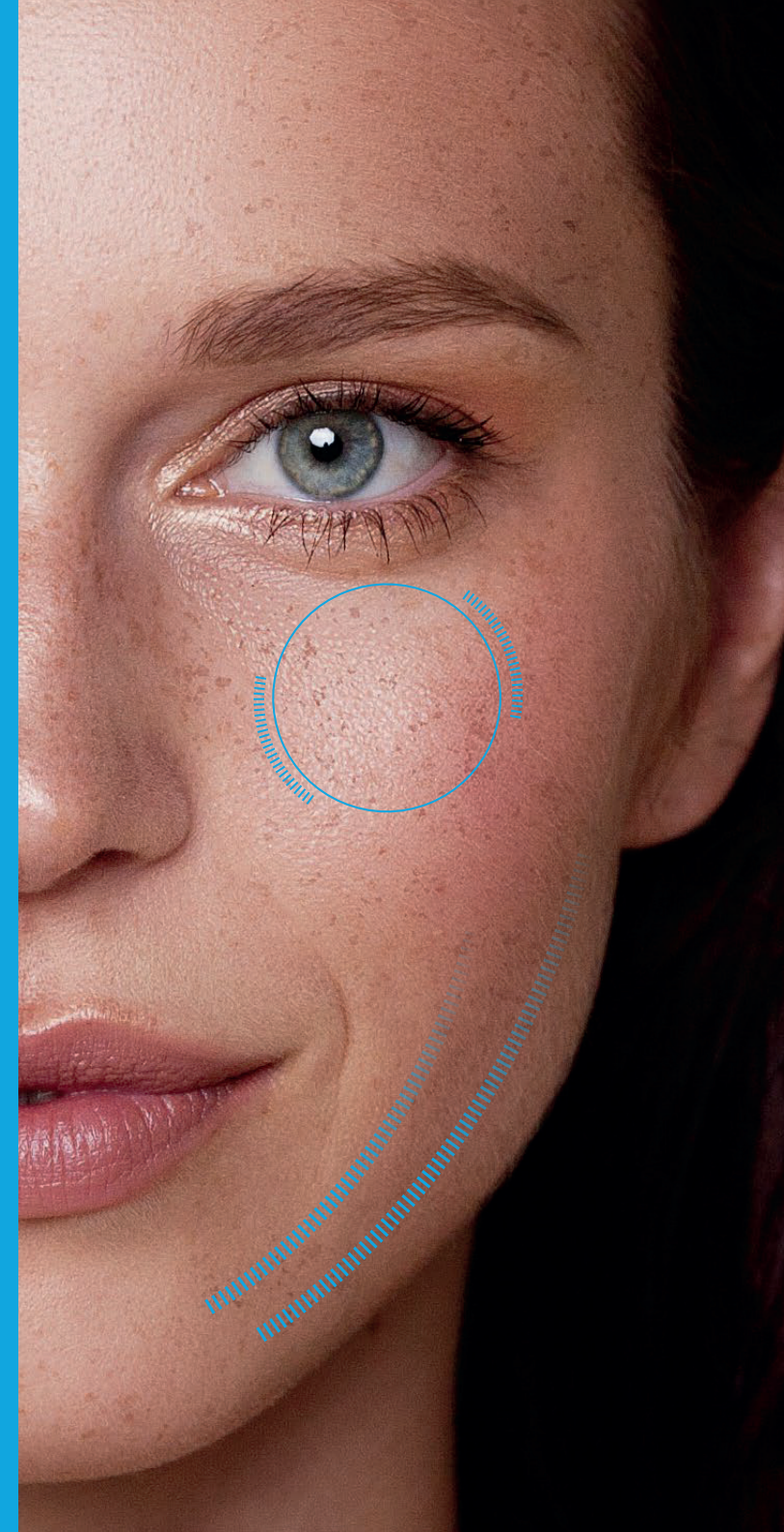

**ENDOLIFT®** IS A MINIMALLY INVASIVE MEDICAL LASER TREATMENT THAT USES EUFOTON® INNOVATIVE LASEMAR® 1500 TO STIMULATE BOTH DEEP AND SUPERFICIAL LAYERS OF THE SKIN, TIGHTEN AND RETRACT THE CONNECTIVE SEPTUM, AND STIMULATE NEW DERMAL COLLAGEN FORMATION.








HOW IT WORKS

The office-based Endolift® treatment requires Eufoton® specific FTF® - Fiber to Fiber system - micro optical fibers, that are easily inserted, without any incisions, under the skin directly in the superfcial hypodermis, creating a micro-tunnel orientated along the anti-gravitational vectors and, after the treatment, the fibers are removed. Eufoton® offers several calibers depending on the area to treat: 200 - 300 - 400 - 600 micron.

While passing through the dermis, these FTF® micro optical fibers act like an intradermal light path and transmit laser energy, offering significant, visible results. The procedure involves minimal to no downtime and it does not have the pain or recovery time that is associated with surgical procedures. Patients can return to work and normal activity within a few hours.

The wavelength of **LASEmaR® 1500** (1470nm) has an ideal interaction with water and fat, which activates neo-collagenesis and the metabolic functions in the extracellular matrix. This results in skin retraction and tightening.

ENDOLIFT® ADVANTAGES

-  Office based procedure.
-  Only one treatment.
-  Safe and immediate results.
-  Long term effect.
-  No incisions.
-  Minimal or no post treatment recovery time.
-  Possibility of therapeutic combination with many surgical and aesthetic treatments.

Results are both immediate and long-term.

The area will continue to improve for several months following the **Endolift®** procedure as additional collagen builds in the deep layers of the skin.


FTF®- FIBER-TO-FIBER SYSTEM

**Endolift®** Micro optical fibers are available in two kinds: bare and radial.



FACE RESULTS


LOWER EYELID before/after

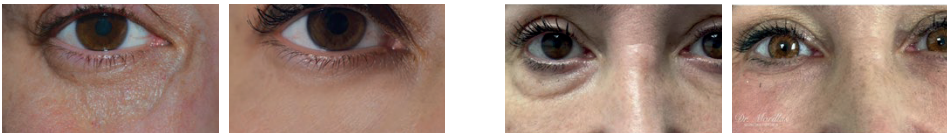

MANDIBULAR CONTOUR - JAWLINE before/after




UNDER CHIN before/after

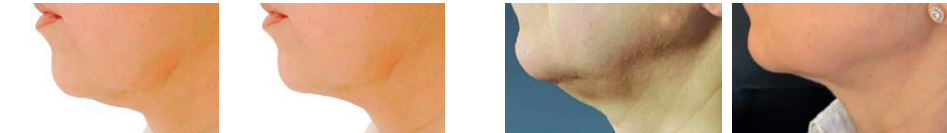

LOWER MIDDLE THIRD OF THE FACE before/after



NECK before/after

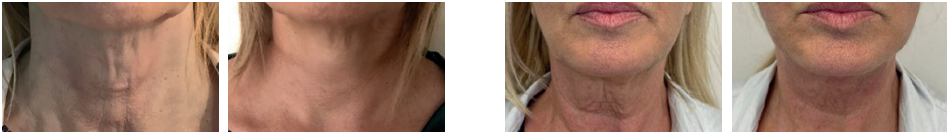

Before and after photographs are courtesy of: M. Baldessin MD, R. Dell'Avanzato MD, A. Dini MD, R. Forte MD, L. Morillas, D. Romano MD, P. Verma MD.