



JOVENA[®]
TOTAL SKIN REGENERATION


brera[®]
medical technologies



JOVENA[®]



TOTAL SKIN REGENERATION

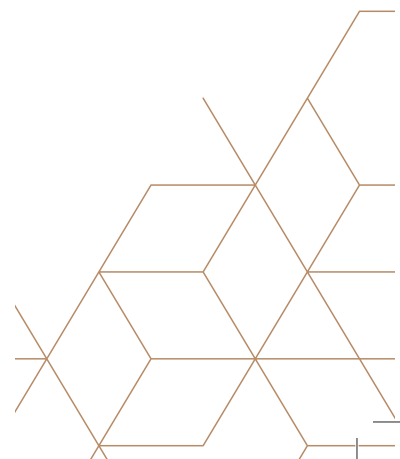



AN INTEGRATIVE
SYSTEM BASED
ON TWO ENERGIES

- ~ RF DIATERMOCONTRACTION®
- ~ FRACTIONAL PLASMA®

FOR SKIN REGENERATION

JOVENA®



BIO-REGENERATIVE PLATFORM

Jovena® takes on a unique approach to skin regeneration by combining two patented modalities in a single device.

Fractional Plasma® for minimally invasive microsurgery and skin resurfacing
Facestim® powered by RF Diatermo Contraction® to reach deeper dermal layers and muscles.

Maximize the activity and profitability of your clinic all year long with unmatched protocols.

Jovena® offers breakthrough treatment options based on your patient's needs.





JOVENA®

TOTAL SKIN REGENERATION




EPIDERMAL
RESURFACING

DEEP DERMAL
REMODELING

MUSCLES
RESTORATION



4 DIFFERENT
MODULES,
HIGHLY
INNOVATIVE
FOR AMAZING
EFFECTS

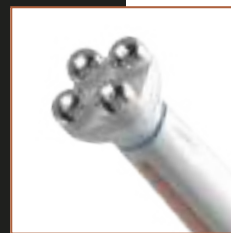

RF bQUAD®



plasmaROLL®




plasmaTIP®



plasmaPRO™

Each **HANDPIECE**
is specifically designed
to maximize treatment results
while minimizing downtime.

JOVENA®
is able to meet your
patients' desired outcomes
by providing comfortable
treatments with little to no pain.

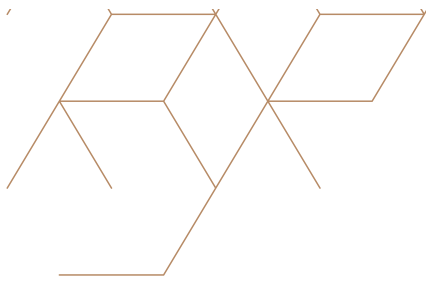

MULTILAYERED TREATMENT OPTIONS FOR A FULL SKIN REGENERATION

- ~ FACE RESURFACING
- ~ DÉCOLLETAGE REVITALIZATION
- ~ MUSCLE TONING
- ~ FACESTIM®
- ~ MINIMALLY INVASIVE MICROSURGERY



JOVENA®






FaceStim[®]

ELEVATE YOUR BEAUTY ROUTINE
Experience the power of facial muscles stimulation

STIMULATING TREATMENT OF INNER MUSCLES

FACESTIM[®] is one of the unique features of Jovena[®]. An intense facial workout helping to tight, volumize, and tone the skin and muscles.

FACESTIM[®] gives the face a healthier and naturally younger look in a non-invasive and comfortable way.

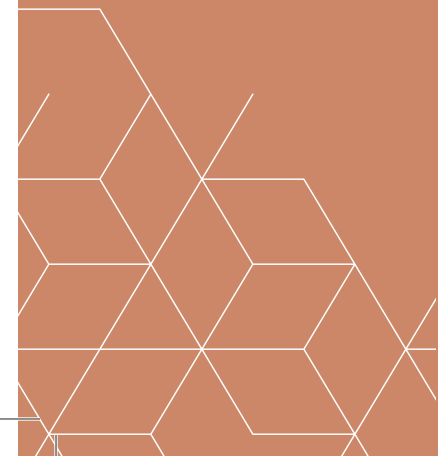


MUSCLES RESTORATION

FaceStim® by Jovena® is based on the simultaneous combination of RF diathermy associated with muscle contractions.

The RF Diathermic Pulse, thanks to its own bio-recognized characteristics, produces massive fibers recruitment while preserving the outstanding comfort during the treatment session.

The muscular contractions generated by the proprietary RF Diathermic Pulse can be combined with the thermic action produced by the same signal to promote skin tightening, firming and toning.

bQUAD® applicator has a patented reticular geometry for optimized and homogenous heat diffusion, a great tool to perform treatments all year long.

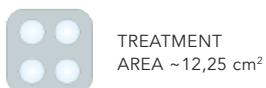

TREATMENT
AREA ~12,25 cm²



2 BIPOLAR
HEMISPHERICAL
ELECTRODES



NO RETURN
ELECTRODE



OPTIMIZED
HEAT
DIFFUSION



Patent n.: 202018000003566


THE EVOLUTION
OF BRERA MEDICAL'S
FRACTIONAL PLASMA®

3 DIFFERENT MODULES FOR
TAILOR-MADE TREATMENTS

plasmaROLL®

plasmaTIP®

plasmaPRO™

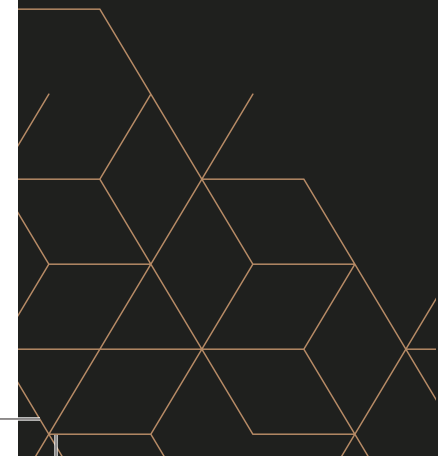

plasmaROLL®

EPIDERMAL RESURFACING

plasmaROLL® is a gentle and safe Fractional Plasma® applicator that targets the superficial layers of the skin and promotes micro-retraction with minimal downtime.


Designed for treatments on large areas of the body, face, neck and décolletage, plasmaROLL® is made of a disposable roll with a dense matrix of micro-holes to assure a fixed distance from the skin.

This also diminishes the risks of undesirable side effects on particularly delicate areas.


plasmaROLL® enables superficial treatments in the epidermis, optimizing epidermal cell renewal.

DAC[™]
D E E P
A C T I O N
C O N T R O L

MICRO-RETRACTION

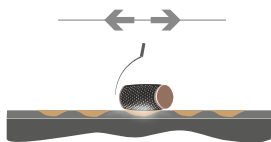

RETRACTION

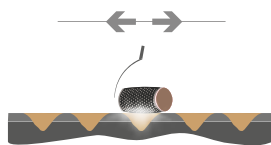

Patent n.: 20201800002599

JOVENA®

plasmaTIP®

SAFE AND CONTROLLED PLASMA DELIVERY

plasmaTIP® provides a precise delivery of Fractional Plasma® to treat common skin aging issues such as active acne, scars and wrinkles.


It is an ideal tool for patients who want to get rid of the aging signs with minimal downtime.

The patented plasmaTIP® is an encapsulated disposable needle tip designed for maximum control of the shot. Its glass acts as a spacer system to keep a fixed distance from the skin and as air booster to assure constant ionization during the session.


The DAC™ (Deep Action Control) provides two modalities: spot and linear techniques. The operator can choose to stay on one fixed area during the shot with the spot option. Selecting the linear option allows the operator a fast sliding treatment with the benefits of automatic shooting.

DAC™
D E E P
A C T I O N
C O N T R O L

SPOT

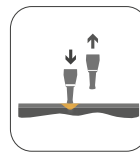

LINEAR

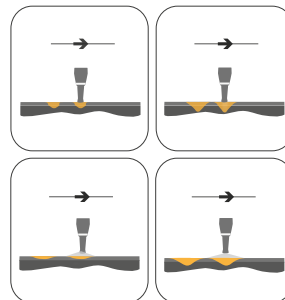

SPRAY

Patent n.: 202019000001575

JOVENA®

MINIMALLY INVASIVE MICROSURGERY

plasmaPRO™ is the latest evolution of leading Plasmage®. It has been enhanced by a selective applicator designed to deliver Fractional Plasma® energy.


The technology creates controlled damages from the epidermis to the superficial dermis with tailored spots in non-ablative or ablative outputs.

The DAC™ treatment modality in the plasmaPRO™ allows the operator to choose among the options super thin, medium or medium-wide spots with deep action and the innovative spray program to target medium size areas with limited damages but with high sublimation.

plasmaPRO™ by Jovena® will continue to be a reference tool for periorbital treatments (blefaroplasma®) and benign lesions by introducing innovative features.

DAC™

D E E P
A C T I O N
C O N T R O L

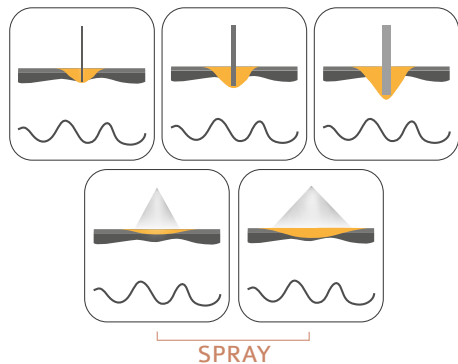

BEFORE /
AFTER



FaceStim[®]

DR. DINKO KALITERNA
DERMATOLOGIST - CROATIA



DR. DANIELA TAHER -
DERMATOLOGIST - ROMANIA

BEFORE /
AFTER

FRACTIONAL
PLASMA®



DR. DINKO KALITERNA -
DERMATOLOGIST - CROATIA



DR. DINKO KALITERNA -
DERMATOLOGIST - CROATIA



BRERA MEDICAL TECHNOLOGIES.
A SUCCESSFUL STORY, ALL MADE IN ITALY.

FOR 40 YEARS, BRERA FAMILY HAS BEEN DEVELOPING AND PRODUCING MEDICAL DEVICES
FOR AESTHETICS, DERMATOLOGY AND PHYSIOTHERAPY WITH A PROFOUNDLY ITALIAN SPIRIT
IN ITS ABILITY TO CREATE, INNOVATE, EXPERIMENT AND INSPIRE A CONSTANTLY EVOLVING MARKET.

brera[®]
medical technologies

VIA ANTONIO DI GENIO 26
84061 OGLIASTRO CILENTO (SA) - ITALY
+39 0974 49.36.516
WWW.BRERAMEDICAL.COM - CONTACT@BRERAMEDICAL.COM

All the trademarks included [®]™ are registered and may not be reproduced.
Jovena[®] device is a registered and has a protected design.
The data, images and all others described here have no contract value, can be modified at any time without prior notice.

